




MS35B

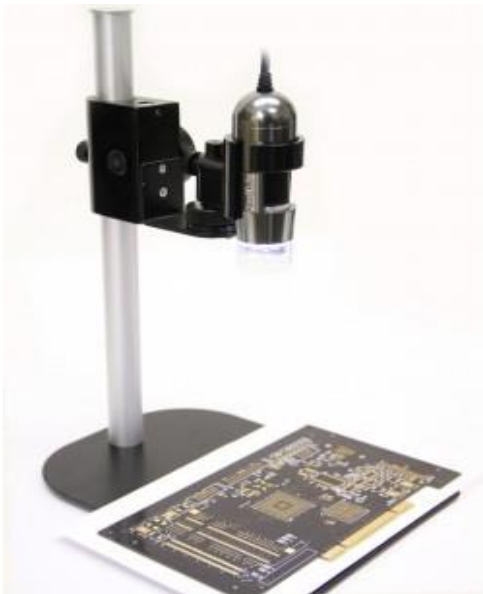

SPECIFICATIONS

Warranty Period	1 year
Manufactured in	China
Service & Support	U.S. (www.dinolite.us) English & Spanish
Regulatory Approval	ROHS
Material	Stand: Metal, Holster: Plastic, Base: Metal
Flexibility	Multi Angle
Dimensions	20.5cm x 12.5cm x 30.5cm (8.07" x 4.92" x 12")
Compatible with	All Dino-Lite Handheld Microscopes
Package Includes	Stand, Holster
Package Dimensions	30.5cm x 15.2cm x 6.7cm (4.92" x 8.07" x 12")
Package Weight	1.2 kg (2.65 lb)
Product Instructions	No
Software Instructions	No

The MS35B is a versatile, lightweight and economical vertical stand with fine adjustment designed for use with the Dino-Lite handheld digital microscope cameras. It has a range of up to 10 inches of vertical working distance from object (up to 20 inches with the optional extension post: MSA2P). The main body assembly with built in fine adjustment slides up or down into position on the vertical post and can be rotated 360 degrees, then secured into position with a rear set screw.

The Dino-Lite handheld microscope (not included) slides directly into a holster which is mounted onto a peg located on the main body. A ball and socket type peg allows for up to 45 degrees of adjustable lateral viewing angles while an optional, separate rigid vertical peg provides a dedicated straight-on view allowing for more accurate measurements. The precision machine crafted components provide smooth adjustment and positioning in a well balanced and stable ergonomic design. Manufactured from a combination of resilient stainless steel and lightweight recycled aluminum, the MS35B can easily and conveniently be moved from one work area to another and will endure the harshest of work environments